

## AKAP10 Polyclonal Antibody

Catalog number: 12356-1-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

AKAP10; D-AKAP2; MGC9414;

PRKA10

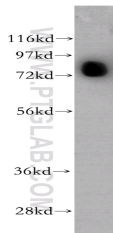

HEK-293 cells were subjected to SDS PAGE followed by western blot with 12356-1-AP (AKAP10 Antibody) at dilution of 1:1500

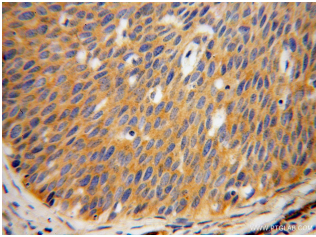

Immunohistochemistry of paraffin-embedded human cervical cancer using 12356-1-AP (AKAP10 Antibody) at Dilution 1:100 (under 40x lens)

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	HEK-293 cells, HeLa cells
Calculated AKAP10 MW:	662aa, 74kd
Observed AKAP10 MW:	74kd
Positive IHC detected in:	Human cervical cancer
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag3020
GenBank accession number:	BC017055
Gene ID (NCBI):	11216
Full name:	A kinase (PRKA) anchor protein 10

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.