

STXBP6 Polyclonal Antibody

Catalog number: 10976-4-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

STXBP6; FLJ39638; HSPC156;

amisyn

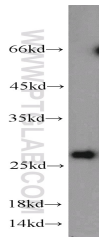

mouse brain tissue were subjected to SDS PAGE followed by western blot with 10976-4-AP(STXBP6 Antibody) at dilution of 1:500

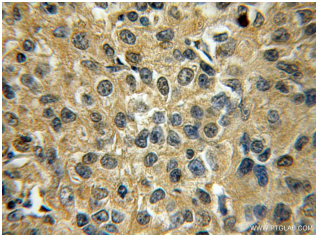

Immunohistochemistry of paraffin-embedded human lung cancer using 10976-4-AP(STXBP6 Antibody) at Dilution 1:50 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Mouse brain tissue, human heart tissue
Calculated STXBP6 MW:	25kd or 45kd
Observed STXBP6 MW:	28-30kd
Positive IHC detected in:	Human lung cancer
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag1427
GenBank accession number:	BC009499
Gene ID (NCBI):	29091
Full name:	Syntaxin binding protein 6 (amisyn)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.