

ACAA2 Polyclonal Antibody

Catalog number: 11111-1-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

ACAA2; DSAEC; FLJ35992;

FLJ95265

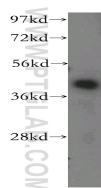

HeLa cells were subjected to SDS PAGE followed by western blot with 11111-1-AP (ACAA2 Antibody) at dilution of 1:500

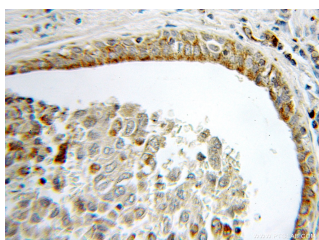

Immunohistochemistry of paraffin-embedded human lung cancer using 11111-1-AP (ACAA2 Antibody) at Dilution 1:100 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Cited applications:	IHC, WB
Species specificity:	Human, mouse, rat; other species not tested.
Cited species:	Rat
Positive WB detected in:	HeLa cells, HEK-293 cells
Calculated ACAA2 MW:	42kd
Observed ACAA2 MW:	42kd
Positive IHC detected in:	Human lung cancer
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP =

Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag1609
GenBank accession number:	BC001918
Gene ID (NCBI):	10449
Full name:	Acetyl-Coenzyme A acyltransferase 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.