

ASB8 Polyclonal Antibody

Catalog number: 11735-1-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

ASB8; FLJ21255; MGC5540;

PP14212

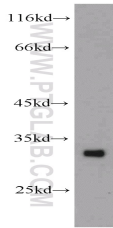

human heart tissue were subjected to SDS PAGE followed by western blot with 11735-1-AP(ASB8 Antibody) at dilution of 1:300

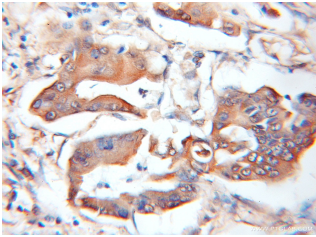

Immunohistochemistry of paraffin-embedded human colon cancer using 11735-1-AP(ASB8 Antibody) at Dilution 1:100 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human heart tissue
Calculated ASB8 MW:	288aa,32kd
Observed ASB8 MW:	32kd
Positive IHC detected in:	Human colon cancer
Recommended dilution:	WB: 1:200-1:1000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag2361
GenBank accession number:	BC001321
Gene ID (NCBI):	140461
Full name:	Ankyrin repeat and SOCS box-containing 8

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.