

PPP5C Polyclonal Antibody

Catalog number: 11715-1-AP

Size: 75 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


PPP5C; FLJ36922, PP T, PP5,

PPP5, PPP5C, PPT, Protein

phosphatase T

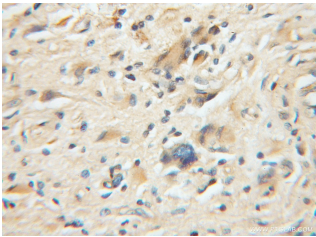

Immunohistochemical of paraffin-embedded human gliomas using 11715-1-AP (PPP5C antibody) at dilution of 1:100 (under 10x lens)

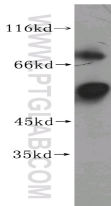

mouse testis tissue were subjected to SDS PAGE followed by western blot with 11715-1-AP (PPP5C antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated PPP5C MW:	499aa, 57 kDa
Observed PPP5C MW:	57 kDa
Positive WB detected in	Mouse testis tissue, human liver tissue
Positive IHC detected in	Human gliomas tissue
Recommended dilution:	WB: 1:1000-1:10000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2288
GenBank accession number:	BC001970
Gene ID (NCBI):	5536
Full name:	Protein phosphatase 5, catalytic subunit

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.