

TROAP Polyclonal Antibody

Catalog number: 13634-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


TROAP; TASTIN, TROAP,

Trophinin assisting protein,

Trophinin associated protein

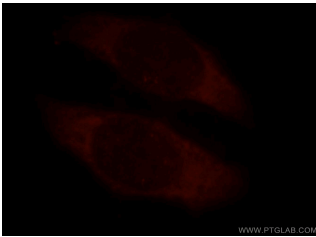

Immunofluorescent analysis of MCF-7 cells, using TROAP antibody 13634-1-AP at 1:25 dilution and Rhodamine-labeled goat anti-rabbit IgG (red).

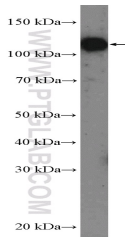

HeLa cells were subjected to SDS PAGE followed by western blot with 13634-1-AP (TROAP Antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, IF, WB
Species specificity:	Human; other species not tested.
Calculated TROAP MW:	778aa,84 kDa
Observed TROAP MW:	110 kDa
Positive WB detected in:	HeLa cells, PC-3 cells
Positive IF detected in:	MCF-7 cells
Recommended dilution:	WB: 1:200-1:1000 IF: 1:10-1:100

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4562
GenBank accession number:	BC032583
Gene ID (NCBI):	10024
Full name:	Trophinin associated protein (tastin)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.