

NECAB2 Polyclonal Antibody

Catalog number: 12257-1-AP

Size: 47 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

NECAB2; EFCBP2, NECAB2,


Stip 2

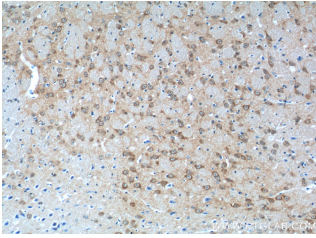

Immunohistochemistry of paraffin-embedded mouse brain tissue slide using 12257-1-AP (NECAB2 Antibody) at dilution of 1:50 (under 10x lens)

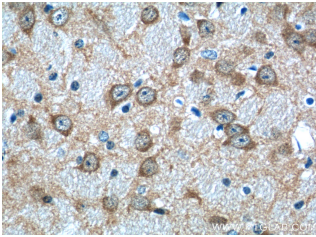

Immunohistochemistry of paraffin-embedded mouse brain tissue slide using 12257-1-AP (NECAB2 Antibody) at dilution of 1:50 (under 40x lens)

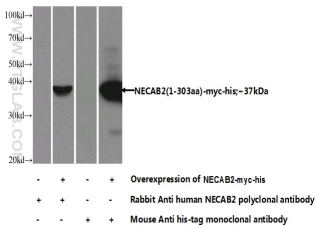

Transfected HEK-293 cells were subjected to SDS PAGE followed by western blot with 12257-1-AP (NECAB2 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated NECAB2 MW:	386aa, 42 kDa
Observed NECAB2 MW:	
Positive WB detected in	Transfected HEK-293 cells
Positive IHC detected in	Mouse brain tissue
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2909
GenBank accession number:	BC016979
Gene ID (NCBI):	54550
Full name:	N-terminal EF-hand calcium binding protein 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.