

## ACYP2 Polyclonal Antibody

Catalog number: 12302-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


ACYP2; ACYP2,

Acylphosphatase 2, ACYP,

ACYM, acylphosphatase 2,  
muscle type

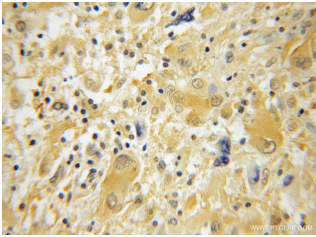

Immunohistochemical of paraffin-embedded human medulloblastoma using 12302-1-AP(ACYP2 antibody) at dilution of 1:50 (under 10x lens)

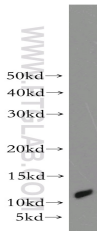

mouse skeletal muscle tissue were subjected to SDS PAGE followed by western blot with 12302-1-AP(ACYP2 antibody) at dilution of 1:500

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Calculated ACYP2 MW:	99aa,11 kDa
Observed ACYP2 MW:	11kd
Positive WB detected in	Mouse skeletal muscle tissue
Positive IHC detected in	Human medulloblastoma
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag2948
GenBank accession number:	BC012290
Gene ID (NCBI):	98
Full name:	Acylphosphatase 2, muscle type

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.