

HMGB4 Polyclonal Antibody

Catalog number: 12787-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

HMGB4; HMGB4, high mobility

group box 4, dJ11007G16.5,

FLJ40388

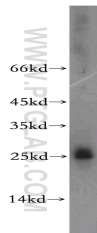

human testis tissue were subjected to SDS PAGE followed by western blot with 12787-1-AP(HMGB4 antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Calculated HMGB4 MW:	186aa, 22 kDa
Observed HMGB4 MW:	22-25kd
Positive WB detected in	Human testis tissue
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3503
GenBank accession number:	BC021180
Gene ID (NCBI):	127540
Full name:	High-mobility group box 4

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.