

PEX5L Polyclonal Antibody

Catalog number: 13084-1-AP


Size: 29 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


PEX5L; Peroxin 5 related protein, PEX2 related protein, PEX5 like protein, PEX5 related protein, PEX5L, PEX5R, Pex5Rp, PXR2, PXR2B

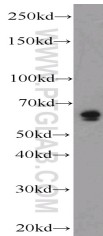

mouse brain tissue were subjected to SDS PAGE followed by western blot with 13084-1-AP(PEX5L antibody) at dilution of 1:600

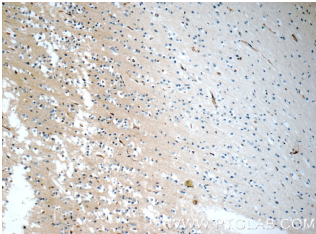

Immunohistochemistry of paraffin-embedded human brain tissue slide using 13084-1-AP (PEX5L Antibody) at dilution of 1:50 (under 10x lens)

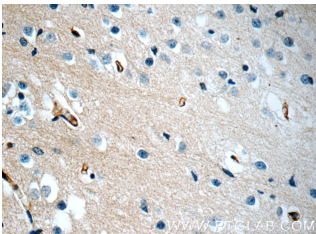

Immunohistochemistry of paraffin-embedded human brain tissue slide using 13084-1-AP (PEX5L Antibody) at dilution of 1:50 (under 40x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated PEX5L MW:	626aa, 70kd
Observed PEX5L MW:	65-70 kDa
Positive WB detected in:	Mouse brain tissue, rat brain tissue
Positive IHC detected in:	Human brain tissue
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3749
GenBank accession number:	BC036183
Gene ID (NCBI):	51555
Full name:	Peroxisomal biogenesis factor 5-like

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.