

## PCDHA6 Polyclonal Antibody

Catalog number: 12853-1-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

PCDHA6; CNR2; CNRN2;

CNRS2; CRNR2; PCDH-ALPHA6

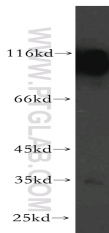

human brain tissue were subjected to SDS PAGE followed by western blot with 12853-1-AP(PCDHA6 Antibody) at dilution of 1:1000

### Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human brain tissue, K-562 cells
Calculated PCDHA6 MW:	103kd
Observed PCDHA6 MW:	75kd, 103kd
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag3493
GenBank accession number:	BC036674
Gene ID (NCBI):	56142
Full name:	Protocadherin alpha 6

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.