

# Futura

+4°C



LABH-30-DTX

## Explosion Proof (Hazardous Locations) Dual Temperature Refrigerator/Freezer

Our **FUTURA Silver Series** 30 cubic foot explosion proof (hazardous locations) laboratory refrigerator/freezer is designed to ensure safety for the storage of volatile materials meeting both NFPA and OSHA guidelines.

### Standard Product Features:

- Manual defrost
- Large 30 cu. ft. interior
- Refrigerator factory preset to 4°C; freezer factory preset to -20°C
- Meets NFPA and OSHA guidelines for Hazardous Locations, Class I, Division II, Group C & D.\*
- Adjustable mechanical thermostat sealed in UL listed explosion proof box
- Must be hard wired in conduit to 120 VAC power source. Power cords with plugs are **not acceptable**.
- Magnetic door gasket for positive seal
- High density urethane foam cabinet and door insulation
- CFC-free refrigerant (R134a), foam insulation and packaging
- Leveling Legs
- White textured cabinet and door
- Heavy duty wire shelves
- Key lock door, right hand swing
- One interior and exterior solid door. Heavy duty, stainless steel pull handle on exterior door
- 115V operation
- UL listed
- Warranties: 2 years parts and labor; 5 years compressor parts only

*\* For questions regarding the correct placement of this unit in your laboratory, please contact your LabRepCo cold storage specialist.*


### Optional Features:

- Casters
- Additional full shelf

# LABRepCo

One Company • Many Solutions


Top View

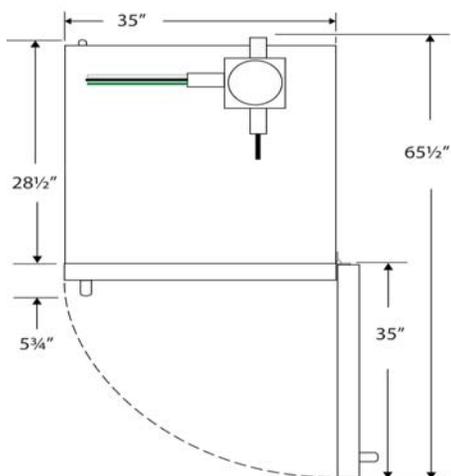

Front Exterior View



Side Exterior View

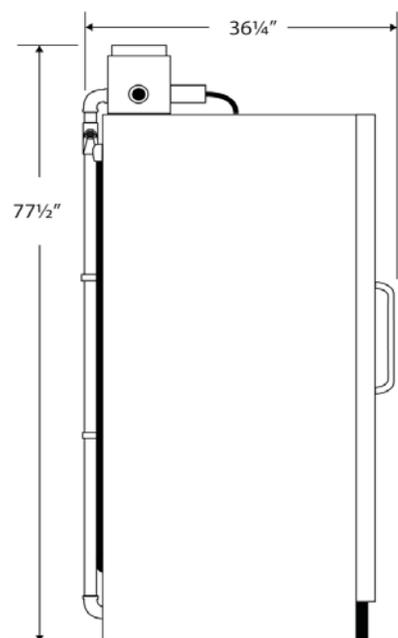

Models	LABH-30-DTX
Weight lbs.	368
Overall Height in.	77 1/2*
Overall Width in.	35 3/4
Overall Depth in.	36 1/4
Interior Refrigerator	
Height in.	50 1/4
Width in.	32
Depth in.	22 1/2
Interior Freezer	
Height in.	12
Width in.	29
Depth in.	22 1/4
Interior Cu. Ft.	
Refrigerator	25
Freezer	5
Number of Shelves	3 adjustable (refrig)
Defrost Type	Manual
Compressor Size HP	1/3
Plug (receptacle)	N/A**
Total Amp Draw	6
Voltage	115V, 1PH, 60Hz

\*Leveling legs are fully recessed. \*\*Must be hard wired in conduit to 120 VAC power source. Power cords with plugs are not acceptable.

Specifications subject to change without notice.



LABH-30-DTX interior

