

gamma catenin Polyclonal Antibody

Catalog number: 11146-1-AP


Size: 40 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

JUP; ARVD12, Catenin gamma, CTNNG, Desmoplakin 3, Desmoplakin III, DP3, DP111, junction plakoglobin, JUP, gamma Catenin, PDGB, PKGB

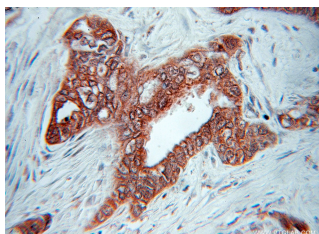

Immunohistochemical of paraffin-embedded human pancreas cancer using 11146-1-AP (JUP antibody) at dilution of 1:50 (under 10x lens)

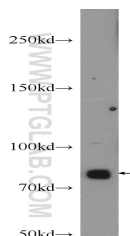

mouse heart tissue were subjected to SDS PAGE followed by western blot with 11146-1-AP (JUP Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, IHC, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated gamma catenin MW:	82kd
Observed gamma catenin MW:	82 kDa
Positive WB detected in	Mouse heart tissue, A431 cells, HeLa cells, human heart tissue
Positive IHC detected in	Human pancreas cancer tissue, human breast cancer tissue, human colon cancer tissue, human skin tissue
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1627
GenBank accession number:	BC011865
Gene ID (NCBI):	3728
Full name:	Junction plakoglobin

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.