

NGB Polyclonal Antibody

Catalog number: 13499-1-AP


Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

NGB;

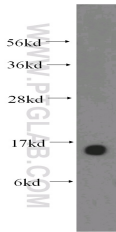

human brain tissue were subjected to SDS PAGE followed by western blot with 13499-1-AP(NGB Antibody) at dilution of 1:1000

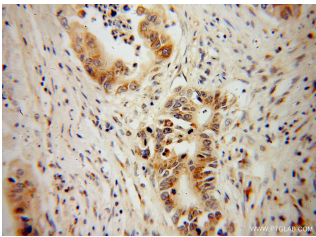

Immunohistochemistry of paraffin-embedded human pancreas cancer using 13499-1-AP(NGB Antibody) at Dilution 1:100 (under 40x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human brain tissue
Calculated NGB MW:	151aa,16kd
Observed NGB MW:	16kd
Positive IHC detected in:	Human pancreas cancer
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag4325
GenBank accession number:	BC032509
Gene ID (NCBI):	58157
Full name:	Neuroglobin

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.