

FCN1 Polyclonal Antibody

Catalog number: 11775-1-AP

Size: 25 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

FCN1; FCN1, FCNM, Ficolin 1,

Ficolin A, Ficolin alpha, M

ficolin

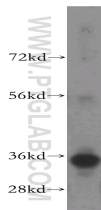

human liver tissue were subjected to SDS PAGE followed by western blot with 11775-1-AP(FCN1 antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB
Species specificity:	Human; other species not tested.
Calculated FCN1 MW:	326aa,35 kDa
Observed FCN1 MW:	35 kDa
Positive WB detected in	Human liver tissue, K-562 cells
Recommended dilution:	WB: 1:1000-1:10000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2444
GenBank accession number:	BC020635
Gene ID (NCBI):	2219
Full name:	Ficolin (collagen/fibrinogen domain containing) 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.