

PMF1 Polyclonal Antibody

Catalog number: 11190-2-AP


Size: 53 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

PMF1;

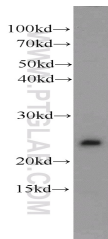

Jurkat cells were subjected to SDS PAGE followed by western blot with 11190-2-AP(PMF1 Antibody) at dilution of 1:500

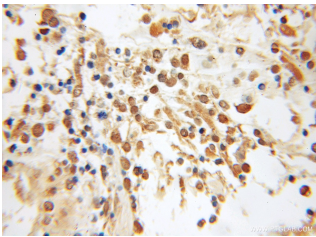

Immunohistochemistry of paraffin-embedded human gliomas using 11190-2-AP(PMF1 Antibody) at Dilution 1:50 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse; other species not tested.
Positive WB detected in:	Jurkat cells, HepG2 cells
Calculated PMF1 MW:	19kd
Observed PMF1 MW:	23kd
Positive IHC detected in:	Human gliomas
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag1672
GenBank accession number:	BC033656
Gene ID (NCBI):	11243
Full name:	Polyamine-modulated factor 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.