

CNTFR Polyclonal Antibody

Catalog number: 10796-1-AP


Size: 23 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

CNTFR; CNTF receptor subunit
alpha, CNTFR, CNTFR alpha

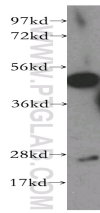

human brain tissue were subjected to SDS PAGE followed by western blot with 10796-1-AP(CNTFR antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated CNTFR MW:	41 kDa
Observed CNTFR MW:	41 kDa
Positive WB detected in	Human brain tissue, mouse cerebellum tissue
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC= Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1189
GenBank accession number:	BC001492
Gene ID (NCBI):	1271
Full name:	Ciliary neurotrophic factor receptor

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.