

## RBP5 Polyclonal Antibody

Catalog number: 11843-1-AP


Size: 50 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


RBP5; CRBP III, CRBP3,  
CRBP111, HRBPiso, RBP5,  
Retinol binding protein 5

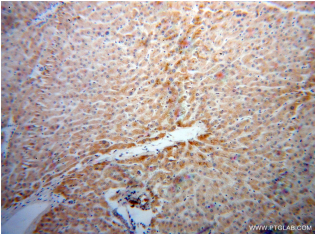

Immunohistochemical of paraffin-embedded human liver cancer using 11843-1-AP(RBP5 antibody) at dilution of 1:50 (under 10x lens)

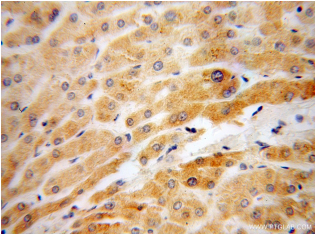

Immunohistochemical of paraffin-embedded human liver cancer using 11843-1-AP(RBP5 antibody) at dilution of 1:50 (under 40x lens)

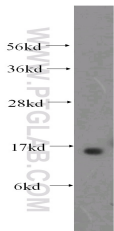

human liver tissue were subjected to SDS PAGE followed by western blot with 11843-1-AP(RBP5 antibody) at dilution of 1:400

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, Pig, Rat; other species not tested.
Calculated RBP5 MW:	135aa,16 kDa
Observed RBP5 MW:	16-18 kDa
Positive WB detected in	Human liver tissue, pig liver tissue, rat liver tissue
Positive IHC detected in	Human liver cancer tissue
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag2485
GenBank accession number:	BC029355
Gene ID (NCBI):	83758
Full name:	Retinol binding protein 5, cellular

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.