

SAFETY DATA SHEET

According to US Regulation 29 CFR 1910.1200 (HazCom 2012)

1. Identification

Product identifier: Sodium Hydroxide, Solution

Other means of identification

Product No.: 0312, 0338, 0895, 3729, 400D, 4690, 5000, 5002, 5661, 5668, 5671, 5672, 5674, 6000, H385, 30400

Recommended restrictions

Recommended use: For Laboratory, Research or Manufacturing Use.

Restrictions on use: Not determined.

Details of the supplier of the safety data sheet

Company Name: Avantor Performance Materials, LLC
Address: 100 Matsonford Rd, Suite 200
Radnor, PA 19087

Telephone: Customer Service: 855-282-6867

Contact Person: Product Information Compliance
E-mail: info@avantormaterials.com

Emergency telephone number:

CHEMTREC: 1-800-424-9300 within US and Canada

2. Hazard(s) identification

Hazard Classification

Physical Hazards

Corrosive to metal Category 1

Health Hazards

Skin Corrosion/Irritation Category 1A

Serious Eye Damage/Eye Irritation Category 1

Specific Target Organ Toxicity -
Single Exposure Category 3¹

Target Organs

1. Respiratory tract irritation.

Unknown toxicity - Health

Acute toxicity, oral	0 %
Acute toxicity, dermal	0 %
Acute toxicity, inhalation, vapor	40 %
Acute toxicity, inhalation, dust or mist	40 %

Environmental Hazards

Acute hazards to the aquatic environment

Category 3

Unknown toxicity - Environment

Acute hazards to the aquatic environment	0 %
Chronic hazards to the aquatic environment	40 %

Label Elements

Hazard Symbol:



Signal Word: Danger

Hazard Statement: May be corrosive to metals.
Causes severe skin burns and eye damage.
May cause respiratory irritation.
Harmful to aquatic life.

Precautionary Statements

Prevention: Keep only in original packaging. Do not breathe dust/mist/vapors. Use only outdoors or in a well-ventilated area. Wear protective gloves/protective clothing/eye protection/face protection. Wash thoroughly after handling. Avoid release to the environment.

Response: Absorb spillage to prevent material damage. IF SWALLOWED: Rinse mouth. Do NOT induce vomiting. IF ON SKIN (or hair): Take off immediately all contaminated clothing. Rinse skin with water [or shower]. Wash contaminated clothing before reuse. IF IN EYES: Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. Immediately call a POISON CENTER/doctor. IF INHALED: Remove person to fresh air and keep comfortable for breathing. Specific treatment (see on this label). Collect spillage.

Storage: Store in a corrosion-resistant container with a resistant inner liner. Store in a well-ventilated place. Keep container tightly closed. Store locked up.

Disposal: Dispose of contents/container to an appropriate treatment and disposal facility in accordance with applicable laws and regulations, and product characteristics at time of disposal.

Hazard(s) not otherwise classified (HNOC): None.

3. Composition/information on ingredients

Mixtures

Chemical Identity	CAS number	Content in percent (%)*
Sodium hydroxide	1310-73-2	20 - 40%

* All concentrations are percent by weight unless ingredient is a gas. Gas concentrations are in percent by volume.

4. First-aid measures

General information:	Get medical advice/attention if you feel unwell. Show this safety data sheet to the doctor in attendance.
Ingestion:	Call a physician or poison control center immediately. Do NOT induce vomiting. If vomiting occurs, keep head low so that stomach content doesn't get into the lungs.
Inhalation:	Move to fresh air. If breathing stops, provide artificial respiration. If breathing is difficult, give oxygen. Call a physician or poison control center immediately.
Skin Contact:	Immediately flush with plenty of water for at least 15 minutes while removing contaminated clothing and shoes. Call a physician or poison control center immediately. Wash contaminated clothing before reuse. Destroy or thoroughly clean contaminated shoes.
Eye contact:	Immediately flush with plenty of water for at least 15 minutes. If easy to do, remove contact lenses. Call a physician or poison control center immediately. In case of irritation from airborne exposure, move to fresh air. Get medical attention immediately.

Most important symptoms/effects, acute and delayed

Symptoms:	Causes severe skin and eye burns. Causes digestive tract burns. Mist or vapor extremely irritating to eyes and respiratory tract.
Hazards:	None known.

Indication of immediate medical attention and special treatment needed

Treatment:	Symptoms may be delayed.
-------------------	--------------------------

5. Fire-fighting measures

General Fire Hazards:	The product is non-combustible. Product is highly caustic. Wear protective gear if spilled during fire fighting.
------------------------------	--

Suitable (and unsuitable) extinguishing media

Suitable extinguishing media:	The product is non-combustible. Use fire-extinguishing media appropriate for surrounding materials.
Unsuitable extinguishing media:	None known.
Specific hazards arising from the chemical:	Product is highly caustic. Wear appropriate protective gear if spilled during firefighting. Contact with metals may evolve flammable hydrogen gas.

Special protective equipment and precautions for firefighters

Special fire fighting procedures:

Move containers from fire area if you can do so without risk. Use water spray to keep fire-exposed containers cool.

Special protective equipment for fire-fighters:

Firefighters must use standard protective equipment including flame retardant coat, helmet with face shield, gloves, rubber boots, and in enclosed spaces, SCBA.

6. Accidental release measures
Personal precautions, protective equipment and emergency procedures:

Put on protective equipment before entering danger area. See Section 8 of the SDS for Personal Protective Equipment. Keep unauthorized personnel away. Keep upwind. Ventilate closed spaces before entering them. Do not touch damaged containers or spilled material unless wearing appropriate protective clothing.

Methods and material for containment and cleaning up:

Neutralize spill area and washings with dilute acetic acid. Absorb spill with vermiculite or other inert material, then place in a container for chemical waste. Collect in a non-combustible container for prompt disposal. Dike far ahead of larger spill for later recovery and disposal.

Notification Procedures:

Inform authorities if large amounts are involved.

Environmental Precautions:

Do not contaminate water sources or sewer. Prevent further leakage or spillage if safe to do so. Avoid discharge into drains, water courses or onto the ground.

7. Handling and storage
Precautions for safe handling:

Use personal protective equipment as required. Avoid breathing mists or vapors. Avoid contact with eyes, skin, and clothing. Do not taste or swallow. Wash hands thoroughly after handling. Do not eat, drink or smoke when using the product. See Section 8 of the SDS for Personal Protective Equipment.

Conditions for safe storage, including any incompatibilities:

Do not store in metal containers. Keep container tightly closed. Store in a well-ventilated place. Store in a dry place.

8. Exposure controls/personal protection
Control Parameters
Occupational Exposure Limits

Chemical Identity	Type	Exposure Limit Values	Source
Sodium hydroxide	Ceiling	2 mg/m3	US. ACGIH Threshold Limit Values (2011)
	Ceil_Time	2 mg/m3	US. NIOSH: Pocket Guide to Chemical Hazards (2010)
	PEL	2 mg/m3	US. OSHA Table Z-1 Limits for Air Contaminants (29 CFR 1910.1000) (02 2006)
	Ceiling	2 mg/m3	US. OSHA Table Z-1-A (29 CFR 1910.1000) (1989)
	Ceiling	2 mg/m3	US. Tennessee. OELs. Occupational Exposure Limits, Table Z1A (06 2008)
	ST ESL	20 µg/m3	US. Texas. Effects Screening Levels (Texas Commission on Environmental Quality) (03 2014)
	AN ESL	2 µg/m3	US. Texas. Effects Screening Levels (Texas Commission on Environmental Quality) (03 2014)
	Ceiling	2 mg/m3	US. California Code of Regulations, Title 8, Section 5155. Airborne Contaminants (08 2010)

Appropriate Engineering Controls No data available.

Individual protection measures, such as personal protective equipment

General information: Good general ventilation (typically 10 air changes per hour) should be used. Ventilation rates should be matched to conditions. If applicable, use process enclosures, local exhaust ventilation, or other engineering controls to maintain airborne levels below recommended exposure limits. If exposure limits have not been established, maintain airborne levels to an acceptable level. An eye wash and safety shower must be available in the immediate work area.

Eye/face protection: Wear safety glasses with side shields (or goggles) and a face shield.

Skin Protection

Hand Protection: Chemical resistant gloves

Other: Wear suitable protective clothing and gloves.

Respiratory Protection: In case of inadequate ventilation use suitable respirator.

Hygiene measures: Provide eyewash station and safety shower. Always observe good personal hygiene measures, such as washing after handling the material and before eating, drinking, and/or smoking. Routinely wash work clothing and protective equipment to remove contaminants. Avoid contact with skin. Do not get in eyes. Wash contaminated clothing before reuse.

9. Physical and chemical properties
--

Appearance

Physical state:	Liquid
Form:	Liquid
Color:	Colorless
Odor:	Odorless
Odor threshold:	No data available.
pH:	14 (20 °C)
Melting point/freezing point:	1 °C
Initial boiling point and boiling range:	115 °C
Flash Point:	Not applicable (water-based)
Evaporation rate:	As water
Flammability (solid, gas):	No data available.
Upper/lower limit on flammability or explosive limits	
Flammability limit - upper (%):	No data available.
Flammability limit - lower (%):	No data available.
Explosive limit - upper (%):	No data available.
Explosive limit - lower (%):	No data available.
Vapor pressure:	As water
Vapor density:	As water
Density:	1.26 g/ml (20 °C)
Relative density:	1.26 (20 °C)
Solubility(ies)	
Solubility in water:	Miscible with water.
Solubility (other):	No data available.
Partition coefficient (n-octanol/water):	No data available.

Auto-ignition temperature: No data available.
Decomposition temperature: No data available.
Viscosity: No data available.

10. Stability and reactivity

Reactivity: Reacts violently with strong acids.

Chemical Stability: Material is stable under normal conditions.

Possibility of hazardous reactions: Hazardous polymerization does not occur.

Conditions to avoid: Avoid contact with oxidizing agents. Reacts violently with strong acids.

Incompatible Materials: Oxidizing agents. Acids. Maleic Anhydride Halogens. Nitromethane. Contact with metals may evolve flammable hydrogen gas.

Hazardous Decomposition Products: Sodium oxides

11. Toxicological information

Information on likely routes of exposure

Inhalation: May cause damage to mucous membranes in nose, throat, lungs and bronchial system.

Skin Contact: Causes severe skin burns.

Eye contact: Causes serious eye damage.

Ingestion: May cause burns of the gastrointestinal tract if swallowed.

Information on toxicological effects

Acute toxicity (list all possible routes of exposure)

Oral Product: No data available.

Dermal Product: No data available.

Inhalation Product: No data available.

Repeated dose toxicity Product: No data available.

Skin Corrosion/Irritation Product: Causes severe skin burns.

Serious Eye Damage/Eye Irritation Product: Causes serious eye damage.

Respiratory or Skin Sensitization

Product: Not a skin sensitizer.

Carcinogenicity

Product: This substance has no evidence of carcinogenic properties.

IARC Monographs on the Evaluation of Carcinogenic Risks to Humans:

No carcinogenic components identified

US. National Toxicology Program (NTP) Report on Carcinogens:

No carcinogenic components identified

US. OSHA Specifically Regulated Substances (29 CFR 1910.1001-1050):

No carcinogenic components identified

Germ Cell Mutagenicity

In vitro

Product: No mutagenic components identified

In vivo

Product: No mutagenic components identified

Reproductive toxicity

Product: No components toxic to reproduction

Specific Target Organ Toxicity - Single Exposure

Product: Respiratory tract irritation.

Specific Target Organ Toxicity - Repeated Exposure

Product: None known.

Target Organs

Specific Target Organ Toxicity - Single Exposure: Respiratory tract irritation.

Aspiration Hazard

Product: Not classified

Other effects: None known.

12. Ecological information

Ecotoxicity:

Acute hazards to the aquatic environment:

Fish

Product: No data available.

Specified substance(s):

Sodium hydroxide
LOAEL (Sander lucioperca, 24 h): ≥ 35 mg/l
LC 50 (Western mosquitofish (Gambusia affinis), 96 h): 125 mg/l
LC 50 (Lepomis macrochirus, 48 h): 99 mg/l

Aquatic Invertebrates

Product: No data available.

Specified substance(s):

Sodium hydroxide
LC 50 (Ophryotrocha diadema, 48 h): 33 - 100 mg/l

LOAEL (Daphnia magna): 40 - 240 mg/l
 LC 50 (Cockle, 48 h): 330 - 1,000 mg/l
 EC 50 (Water flea (Ceriodaphnia dubia), 48 h): 34.59 - 47.13 mg/l
 EC 50 (Ceriodaphnia sp., 48 h): 40.4 mg/l

Chronic hazards to the aquatic environment:

Fish

Product: No data available.

Aquatic Invertebrates

Product: No data available.

Toxicity to Aquatic Plants

Product: No data available.

Persistence and Degradability

Biodegradation

Product: Expected to be readily biodegradable.

BOD/COD Ratio

Product: No data available.

Bioaccumulative potential

Bioconcentration Factor (BCF)

Product: No data available on bioaccumulation.

Partition Coefficient n-octanol / water (log Kow)

Product: No data available.

Mobility in soil:

The product is water soluble and may spread in water systems.

Other adverse effects:

Harmful to aquatic organisms. The product may affect the acidity (pH-factor) in water with risk of harmful effects to aquatic organisms.

13. Disposal considerations

Disposal instructions:

Discharge, treatment, or disposal may be subject to national, state, or local laws. Since emptied containers retain product residue, follow label warnings even after container is emptied.

Contaminated Packaging:

Since emptied containers retain product residue, follow label warnings even after container is emptied.

14. Transport information

DOT

UN Number:	UN 1824
UN Proper Shipping Name:	Sodium hydroxide solution
Transport Hazard Class(es)	
Class:	8
Label(s):	8
Packing Group:	II
Marine Pollutant:	No

Special precautions for user: Keep away from acids.

IMDG

UN Number: UN 1824
 UN Proper Shipping Name: SODIUM HYDROXIDE SOLUTION
 Transport Hazard Class(es)
 Class: 8
 Label(s): 8
 EmS No.: F-A, S-B
 Packing Group: II
 Marine Pollutant: No
 Special precautions for user: Keep away from acids.

IATA

UN Number: UN 1824
 Proper Shipping Name: Sodium hydroxide solution
 Transport Hazard Class(es):
 Class: 8
 Label(s): 8
 Packing Group: II
 Marine Pollutant: No
 Special precautions for user: Keep away from acids.

15. Regulatory information

US Federal Regulations

TSCA Section 12(b) Export Notification (40 CFR 707, Subpt. D)
 None present or none present in regulated quantities.

US. OSHA Specifically Regulated Substances (29 CFR 1910.1001-1050)
 None present or none present in regulated quantities.

CERCLA Hazardous Substance List (40 CFR 302.4):

<u>Chemical Identity</u>	<u>Reportable quantity</u>
Sodium hydroxide	1000 lbs.

Superfund Amendments and Reauthorization Act of 1986 (SARA)

Hazard categories

- Corrosive to metal
- Skin Corrosion or Irritation
- Serious eye damage or eye irritation
- Specific target organ toxicity (single or repeated exposure)

SARA 302 Extremely Hazardous Substance

None present or none present in regulated quantities.

SARA 304 Emergency Release Notification

<u>Chemical Identity</u>	<u>Reportable quantity</u>
Sodium hydroxide	1000 lbs.

SARA 311/312 Hazardous Chemical

<u>Chemical Identity</u>	<u>Threshold Planning Quantity</u>
Sodium hydroxide	10000 lbs.

SARA 313 (TRI Reporting)

None present or none present in regulated quantities.

Clean Air Act (CAA) Section 112(r) Accidental Release Prevention (40 CFR 68.130):

None present or none present in regulated quantities.

Clean Water Act Section 311 Hazardous Substances (40 CFR 117.3):

<u>Chemical Identity</u>	<u>Reportable quantity</u>
Sodium hydroxide	Reportable quantity: 1000 lbs.

US State Regulations

US. California Proposition 65

No ingredient regulated by CA Prop 65 present.

US. New Jersey Worker and Community Right-to-Know Act

<u>Chemical Identity</u>
Sodium hydroxide

US. Massachusetts RTK - Substance List

<u>Chemical Identity</u>
Sodium hydroxide

US. Pennsylvania RTK - Hazardous Substances

<u>Chemical Identity</u>
Sodium hydroxide

US. Rhode Island RTK

<u>Chemical Identity</u>
Sodium hydroxide

International regulations

Montreal protocol

Not applicable

Stockholm convention

Not applicable

Rotterdam convention

Not applicable

Kyoto protocol


Not applicable

Inventory Status:

Australia AICS:	On or in compliance with the inventory
Canada DSL Inventory List:	On or in compliance with the inventory
EINECS, ELINCS or NLP:	On or in compliance with the inventory
Japan (ENCS) List:	On or in compliance with the inventory
Korea Existing Chemicals Inv. (KECI):	On or in compliance with the inventory
Philippines PICCS:	On or in compliance with the inventory
US TSCA Inventory:	On or in compliance with the inventory
New Zealand Inventory of Chemicals:	On or in compliance with the inventory
China Inv. Existing Chemical Substances:	On or in compliance with the inventory
Mexico INSQ:	On or in compliance with the inventory
Taiwan Chemical Substance Inventory:	On or in compliance with the inventory

16. Other information, including date of preparation or last revision

NFPA Hazard ID

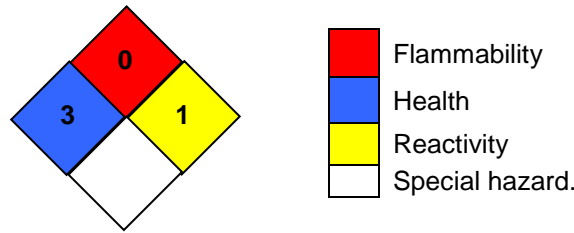

Hazard rating: 0 - Minimal; 1 - Slight; 2 - Moderate; 3 - Serious; 4 - Severe; RNP - Rating not possible

Issue Date:	09-28-2018
Revision Information:	Not relevant.
Version #:	1.3
Source of information:	Sources of information used in preparing this SDS included one or more of the following: results from in house or supplier toxicology studies, information from the Toxicology Data Network (TOXNET), European Chemical Agency (ECHA) substance dossiers, IARC Monographs, US National Toxicology Program data, the Agency for Toxic Substances and Disease Registry, other manufacturer’s SDSs and other sources, as appropriate.
Further Information:	No data available.
Disclaimer:	The information provided in this Safety Data Sheet (SDS) was prepared based on data believed to be accurate as of the date of this SDS. TO THE GREATEST EXTENT PERMITTED BY LAW, AVANTOR PERFORMANCE MATERIALS (“AVANTOR”) EXPRESSLY DISCLAIMS ANY AND ALL REPRESENTATIONS AND WARRANTIES REGARDING THE INFORMATION CONTAINED HEREIN INCLUDING, WITHOUT LIMITATION, AS TO ACCURACY, COMPLETENESS, FITNESS FOR PURPOSE OR USE, MERCHANTABILITY, NON-INFRINGEMENT, PERFORMANCE, SAFETY, SUITABILITY AND STABILITY. This SDS is intended as a guide to the appropriate use, handling, storage and disposal of the product to which it relates by properly trained personnel, and is not intended to be comprehensive. Users of Avantor’s products are advised to perform their own tests and to exercise their own judgment to determine the safety, suitability and appropriate use, handling, storage and disposal of each product and product combination for their own purposes and uses. TO THE GREATEST EXTENT PERMITTED BY LAW, AVANTOR DISCLAIMS LIABILITY FOR, AND BY USING AVANTOR’S PRODUCTS PURCHASER AGREES THAT UNDER NO CIRCUMSTANCES SHALL AVANTOR BE LIABLE FOR, SPECIAL, INDIRECT, INCIDENTAL , PUNITIVE OR CONSEQUENTIAL DAMAGES OF ANY TYPE OR KIND, INCLUDING WITHOUT LIMITATION, FOR LOSS OF PROFITS, REPUTATIONAL DAMAGE, PRODUCT RECALL OR BUSINESS INTERRUPTION.