

DAZL Polyclonal Antibody

Catalog number: 12633-1-AP

Size: 26 µg/150 µl

Source: Rabbit

Isotype: IgG


Synonyms:

DAZL; DAZ homolog, DAZ like
autosomal, DAZH, DAZL,
DAZL1, DAZLA, deleted in
azoospermia like, Deleted in
azoospermia like 1, SPGY like
autosomal, SPGYLA


Applications

Tested applications:	ELISA, WB, IP
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated DAZL MW:	295aa, 33 kDa
Observed DAZL MW:	35kd
Positive WB detected in	Mouse testis tissue
Positive IP detected in	Mouse testis tissue
Recommended dilution:	WB: 1:500-1:5000 IP: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

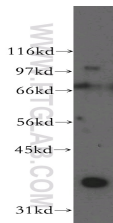

mouse testis tissue were subjected to SDS PAGE followed by western blot with 12633-1-AP (DAZL antibody) at dilution of 1:400

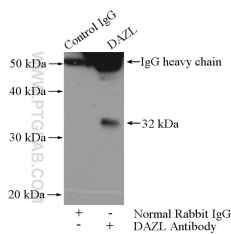

IP Result of anti-DAZL (IP:12633-1-AP, 4µg; Detection:12633-1-AP 1:500) with mouse testis tissue lysate 4400µg.

Immunogen information

Immunogen:	Ag3326
GenBank accession number:	BC027595
Gene ID (NCBI):	1618
Full name:	Deleted in azoospermia-like

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.