

TMOD3 Polyclonal Antibody

Catalog number: 12521-1-AP

Size: 68 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:


TMOD3; TMOD3,

Tropomodulin 3, tropomodulin

3 (ubiquitous), U Tmod,

Ubiquitous tropomodulin,

UTMOD

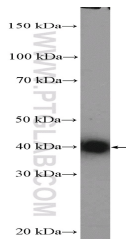

Jurkat cells were subjected to SDS PAGE followed by western blot with 12521-1-AP (TMOD3 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated TMOD3 MW:	352aa, 40 kDa
Observed TMOD3 MW:	40 kDa
Positive WB detected in	Jurkat cells, human heart tissue, human kidney tissue, human lung tissue, mouse brain tissue, mouse kidney tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3215
GenBank accession number:	BC020542
Gene ID (NCBI):	29766
Full name:	Tropomodulin 3 (ubiquitous)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.