

RPN1 Polyclonal Antibody

Catalog number: 12894-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

RPN1; DKFZp686B16177, hrd2,

OST1, RBPH1, Ribophorin 1,

ribophorin I, RPN I, RPN1

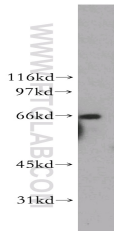

HeLa cells were subjected to SDS PAGE followed by western blot with 12894-1-AP (RPN1 antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB
Cited applications:	WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated RPN1 MW:	69 kDa
Observed RPN1 MW:	66-69 kDa
Positive WB detected in:	HeLa cells, K-562 cells
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3576
GenBank accession number:	BC010839
Gene ID (NCBI):	6184
Full name:	Ribophorin I

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.