

IL32 Polyclonal Antibody

Catalog number: 11079-1-AP


Size: 45 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

IL32; IL 32, IL 32alpha, IL 32beta, IL 32delta, IL 32gamma, IL32, interleukin 32, Natural killer cells protein 4, NK4, TAIF, TAIFa, TAIFb, TAIFc, TAIFd

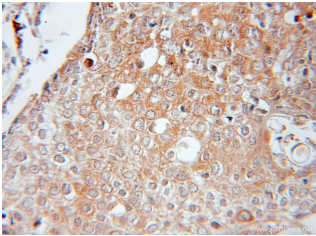

Immunohistochemical of paraffin-embedded human lung cancer using 11079-1-AP(IL32 antibody) at dilution of 1:50 (under 10x lens)

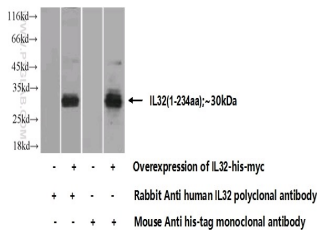

Transfected HEK-293 cells were subjected to SDS PAGE followed by western blot with 11079-1-AP(IL32 Antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, IHC, WB
Species specificity:	Human; other species not tested.
Calculated IL32 MW:	27 kDa
Observed IL32 MW:	
Positive WB detected in	Transfected HEK-293 cells
Positive IHC detected in	Human lung cancer tissue
Recommended dilution:	WB: 1:500-1:5000
	IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1484
GenBank accession number:	BC009401
Gene ID (NCBI):	9235
Full name:	Interleukin 32

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.