

ARHGAP15 Polyclonal Antibody

Catalog number: 11940-1-AP


Size: 22 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

ARHGAP15; BM046

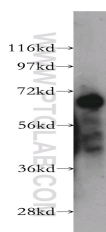

human liver tissue were subjected to SDS PAGE followed by western blot with 11940-1-AP (ARHGAP15 Antibody) at dilution of 1:1000

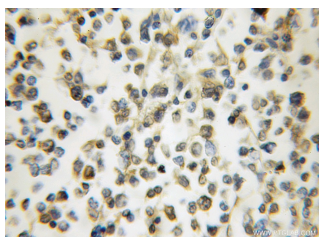

Immunohistochemistry of paraffin-embedded human lymphoma using 11940-1-AP (ARHGAP15 Antibody) at Dilution 1:50 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	Human liver tissue
Calculated ARHGAP15 MW:	475aa, 55kd
Observed ARHGAP15 MW:	65-70kd
Positive IHC detected in:	Human lymphoma
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag2538
GenBank accession number:	BC016701
Gene ID (NCBI):	55843
Full name:	Rho GTPase activating protein 15

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.