

MYF6 Polyclonal Antibody

Catalog number: 11754-1-AP

Size: 45 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

MYF6; bHLHc4, MRF4, Myf 6,

MYF6, Myogenic factor 6,

myogenic factor 6 (herculin)

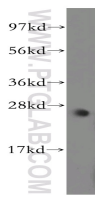

mouse skeletal muscle tissue were subjected to SDS PAGE followed by western blot with 11754-1-AP(MYF6 antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated MYF6 MW:	242aa,27 kDa
Observed MYF6 MW:	27 kDa
Positive WB detected in	Mouse skeletal muscle tissue
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2350
GenBank accession number:	BC017834
Gene ID (NCBI):	4618
Full name:	Myogenic factor 6 (herculin)

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.