

Certified Reference Material

Reference material Certificate

Cannabis Pesticides - AOAC Set 1

Product no.: 56519
Lot no.: BCCG7033
Description of CRM: Pesticides Standard, 5 components in Acetonitrile
Expiry date: MAR 2024
Storage: Store at -20°C ± 5°C
Density (certified) at 20°C: 782.3 kg m⁻³ ± 0.4 kg m⁻³

TraceCERT®
 Traceable Certified Reference Materials

Constituent

Certified values at 20°C and ± Expanded uncertainties, $U = k \cdot u$ ($k = 2$)^{[1][2]}

Constituent	Certified values at 20°C and ± Expanded uncertainties, $U = k \cdot u$ ($k = 2$) ^{[1][2]}	Certified values at 20°C and ± Expanded uncertainties, $U = k \cdot u$ ($k = 2$) ^{[1][2]}
Ancymidol	127.8 mg/kg ± 15.9 mg/kg	100.0 mg/l ± 12.5 mg/l
Chlormequat chloride	129.1 mg/kg ± 3.2 mg/kg	101.0 mg/l ± 2.5 mg/l
Ethephon	126.2 mg/kg ± 7.5 mg/kg	98.7 mg/l ± 5.9 mg/l
Flurprimidol	126.5 mg/kg ± 6.1 mg/kg	99.0 mg/l ± 4.8 mg/l
Phosmet-oxon	127.9 mg/kg ± 6.2 mg/kg	100.1 mg/l ± 4.9 mg/l

Metrological traceability: Traceable to NIST SRM. Details see "Certification process details".^[3]

Measurement method: The certified value is established by weighing of well characterized starting materials in accordance with ISO/IEC 17025^[4]

Intended use: Use this certified reference material (CRM) as a calibrant for pesticide testing in cannabis or any other appropriate matrix.

Instructions for handling and correct use: Shake well and allow to warm to room temperature before opening the ampule. Open a new ampule for each measurement series.

Health and safety information: Please refer to the Safety Data Sheet for detailed information about the nature of any hazard and appropriate precautions to be taken.


Packaging: Amber glass ampule

Accreditation: Sigma-Aldrich Production GmbH is accredited by the Swiss Accreditation Service SAS as reference material producer under no. SRMS 0001 in accordance with international standard ISO 17034^[5].

Certificate issue date: April 26, 2023



H. Sprecher – CRM Operations

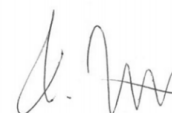

Dr. P. Zell – Approving Officer



ISO 17034
 SRMS 0001



Certification process details:

To guarantee top reliability of the values for this *TraceCERT*[®] certified reference material, well-established procedures were followed. The values have to agree in the range of their uncertainties, but the value from the gravimetric preparation has been chosen as certified value:

1. Gravimetric preparation using well-characterized materials is a practical realization of concentration units, through conversion of mass to amount of substance [6]. If the purity of the materials is demonstrated and if contamination and loss of material is strictly prevented this approach allows highest accuracy. The certified value of *TraceCERT*[®] reference materials is based on this approach and directly traceable to the SI unit kilogram. Therefore, comprehensively characterized materials are used. All balances are certified by DKD and calibrated with OIML Class E2 (up to 12 kg) and F2 (up to 64 kg) weights.
2. The starting material is measured against a certified reference material which is traceable to SRM[®] from NIST followed by gravimetric preparation using balances calibrated with SI-traceable weights. Consequently, the value calculated by this unbroken chain of comparisons is traceable to the reference to which the starting material is compared. Due to the nature of qNMR measurements, different analytes can be characterized using the same SRM[®] from NIST (see table „Details on starting material“ on page 3).
3. For the Chlormequat chloride value, the content of the Chlormequat cation was measured and calculated with a equimolar amount of chloride to the certified content of Chlormequat chloride.
4. The mix in its final packaged form is measured by different LC methods to underpin the gravimetric values.
5. Density measurement is performed in accordance with ISO/IEC 17025^[4]
6. The certificate is designed in accordance with ISO Guide 31. [7]

Homogeneity assessment:

Due to the nature of the production process, a homogeneous solution derives. Nevertheless, homogeneity is assessed by LC measurements and a homogeneity contribution is included into the calculation of content uncertainty of this CRM.


Stability assessment:

A stability study is performed with samples which are stored at different temperatures. The material is tested by LC-MS after certain time intervals and a stability component included in the overall uncertainty.


All sample and reference solutions are mixed with an internal standard solution according to "Guidelines for achieving high accuracy in IDMS"^[8] to give matching signal ratios. When available isotopologues of the corresponding compound are used, for all other compounds closely related isotopologues of the given mixture serve as internal standard. The mixtures are prepared gravimetrically. The analysis of the sample is bracketed by the reference and repeated several times to minimize instrumental variations. The repeatability is included in the combined measurement uncertainty. The set of compounds is analyzed by one LC-MS method. The method was optimized and validated for this set of compounds.

Uncertainty evaluation:

The uncertainty contributions are illustrated by the following cause-effect diagram:

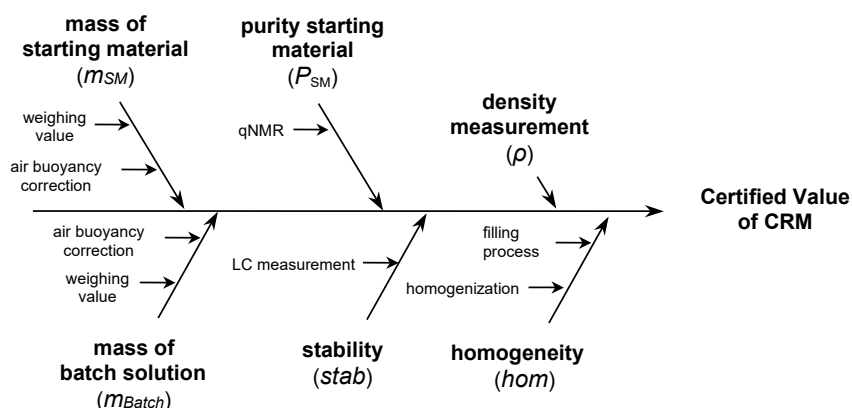

The combined standard uncertainty is calculated by combination of the standard uncertainties of the input estimates according to Eurachem/CITAC Guide "Quantifying Uncertainty in Analytical Measurement" and ISO 17034. [2][5] Expanded uncertainty is then calculated to a confidence level of 95%, typically by multiplying with a coverage factor of $k=2$.

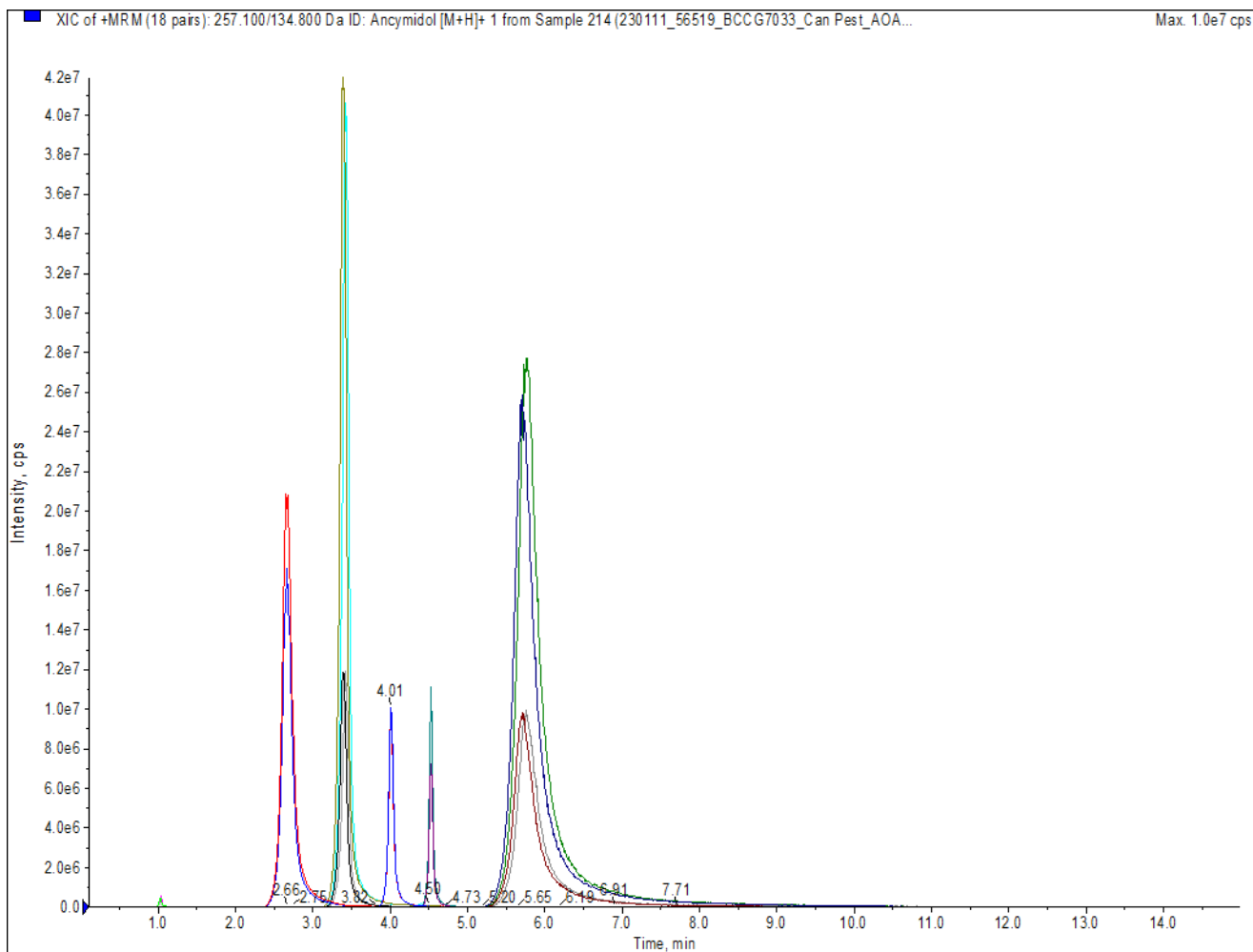

Column: Thermo Trinity P2 (3.0 x 100 mm, 3.0 μm)
 Temperature: 45 °C.
 Mobile Phase: 100 mM ammonium formate in water (A), acetonitrile (B)
 Flow Rate: 0.6 ml/min

Gradient: 0.0 min (95 % A), 1 min (95 % A), 10 min (2 % A), 14.9 min (2 % A), 15 min (95 % A)
 Injection Vol. : 0.2 μL.
 Pre-runtime: 5 min
 *Signal at 2.7 min is internal standard.

Details on starting materials:

#	Compound	CAS Number	t _R / min	Method	NIST SRM traceability
1	Ancymidol	12771-68-5	4.0	qNMR	NIST SRM 84I (KHP)
2	Chlormequat chloride	999-81-5	5.8	qNMR	NIST SRM 84I (KHP)
3	Ethephon	16672-87-0	1.0	qNMR	NIST SRM 84I (KHP)
4	Flurprimidol	56425-91-3	4.5	qNMR	NIST PS1 (Benzoic acid)
5	Phosmet-oxon	3735-33-9	3.4	Mass balance	SI unit kg

References:

- [1] ISO Guide 35:2017, "Reference materials - Guidance for characterization and assessment of homogeneity and stability"
- [2] Eurachem/CITAC Guide, 3rd Ed. (2012), "Quantifying uncertainty in analytical measurement"
- [3] Eurachem/CITAC Guide, 2nd Ed. (2019), "Metrological traceability in chemical measurement" ISO
- [4] The accredited testing laboratory STS 0490 performs the measurements and weighing steps for the certification of this CRM under ISO/IEC 17025:2017, "General requirements for the competence of testing and calibration laboratories"
- [5] 17034:2016, "General requirements for the competence of reference material producers"
- [6] Reichmuth, A., Wunderli, S., Weber, M., Meyer, V. R. (2004), "The uncertainty of weighing data obtained with electronic analytical balances", Microchimica Acta 148: 133-141.
- [7] ISO Guide 31:2015, "Reference materials - Contents of certificates, labels and accompanying documentation"
- [8] Sargent, M., Harrington, C., Harte, R. (eds.), "Guidelines for achieving high accuracy in IDMS", LGC (RSC), Cambridge, 2002.

Certificate of analysis revision history:

Certificate version	Date	Reason for version
01	April 26, 2023	Initial version

Disclaimer:

The purchaser must determine the suitability of this product for its particular use. Sigma-Aldrich Production GmbH makes no warranty of any kind, express or implied, other than its products meet all quality control standards set by Sigma-Aldrich Production GmbH. We do not guarantee that the product can be used for a special application.

