

BCL7B Polyclonal Antibody

Catalog number: 11740-1-AP


Size: 37 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

BCL7B;

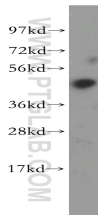

HeLa cells were subjected to SDS PAGE followed by western blot with 11740-1-AP(BCL7B Antibody) at dilution of 1:1000

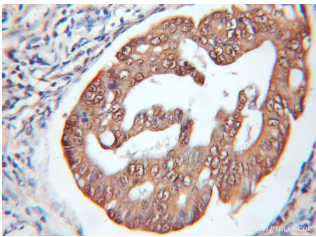

Immunohistochemistry of paraffin-embedded human colon cancer using 11740-1-AP(BCL7B Antibody) at Dilution 1:100 (under 10x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	HeLa cells, human brain tissue
Calculated BCL7B MW:	22kd
Observed BCL7B MW:	45kd
Positive IHC detected in:	Human colon cancer
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag2333
GenBank accession number:	BC000956
Gene ID (NCBI):	9275
Full name:	B-cell CLL/lymphoma 7B

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.