

EEF2K Polyclonal Antibody

Catalog number: 13510-1-AP

Size: 52 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

EEF2K; eEF 2 kinase, eEF 2K,


EEF2K, HSU93850

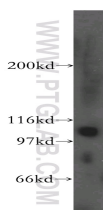

PC-3 cells were subjected to SDS PAGE followed by western blot with 13510-1-AP(EEF2K antibody) at dilution of 1:400

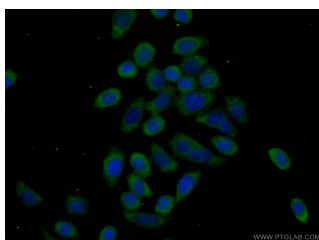

Immunofluorescent analysis of (-20% Ethanol) fixed PC-3 cells using 13510-1-AP(EEF2K Antibody) at dilution of 1:50 and Alexa Fluor 488-conjugated AffiniPure Goat Anti-Rabbit IgG(H+L)

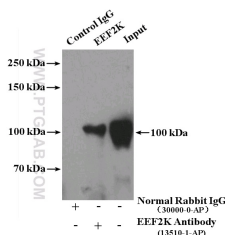

IP Result of anti-EEF2K (IP:13510-1-AP, 4ug; Detection:13510-1-AP 1:500) with PC-3 cells lysate 3200ug.

Applications

Tested applications:	ELISA, WB, IP, IF
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated EEF2K MW:	725aa,100 kDa
Observed EEF2K MW:	100 kDa
Positive WB detected in	PC-3 cells
Positive IP detected in	PC-3 cells
Positive IF detected in	PC-3 cells
Recommended dilution:	WB: 1:500-1:5000 IP: 1:200-1:2000 IF: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4432
GenBank accession number:	BC032665
Gene ID (NCBI):	29904
Full name:	Eukaryotic elongation factor-2 kinase

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.