

CXCL11 Polyclonal Antibody

Catalog number: 10707-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

CXCL11; H174; I-TAC; IP-9; IP9;

MGC102770; SCYB11; SCYB9B;

b-R1

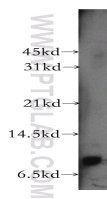

K-562 cells were subjected to SDS PAGE followed by western blot with 10707-1-AP (CXCL11 Antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human; other species not tested.
Positive WB detected in:	K-562 cells
Calculated CXCL11 MW:	10kd
Observed CXCL11 MW:	10 kDa
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag1087
GenBank accession number:	BC005292
Gene ID (NCBI):	6373
Full name:	Chemokine (C-X-C motif) ligand 11

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.