

PPEF1 Polyclonal Antibody

Catalog number: 12875-1-AP

Size: 40 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

PPEF1; PP7, PPEF, PPEF 1,

PPEF1, PPP7C

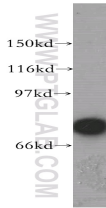

HEK-293 cells were subjected to SDS PAGE followed by western blot with 12875-1-AP(PPEF1 antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated PPEF1 MW:	653aa, 76 kDa
Observed PPEF1 MW:	76 kDa
Positive WB detected in:	HEK-293 cells, human brain tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3553
GenBank accession number:	BC036026
Gene ID (NCBI):	5475
Full name:	Protein phosphatase, EF-hand calcium binding domain 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.