

OSCP1 Polyclonal Antibody

Catalog number: 12598-1-AP

Size: 29 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

C1orf102; C1orf102, hOSCP1,

NOR1, OSCP1, Protein OSCP1

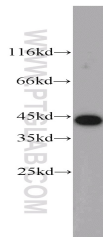

mouse testis tissue were subjected to SDS PAGE followed by western blot with 12598-1-AP(C1orf102 antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated OSCP1 MW:	379aa,43 kDa
Observed OSCP1 MW:	43kd
Positive WB detected in	Mouse testis tissue, mouse brain tissue, mouse lung tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3209
GenBank accession number:	BC018069
Gene ID (NCBI):	127700
Full name:	Chromosome 1 open reading frame 102

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.