

PDSS2 Polyclonal Antibody

Catalog number: 13544-1-AP

Size: 33 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

PDSS2; C6orf210; DLP1;

bA59I9.3; hDLP1

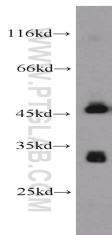

SGC-7901 cells were subjected to SDS PAGE followed by western blot with 13544-1-AP (PDSS2 Antibody) at dilution of 1:600

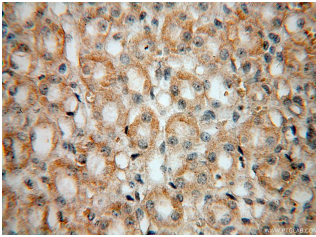

Immunohistochemistry of paraffin-embedded human kidney using 13544-1-AP (PDSS2 Antibody) at Dilution 1:50 (under 40x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	SGC-7901 cells, A375 cells, human testis tissue
Calculated PDSS2 MW:	399aa, 44kd
Observed PDSS2 MW:	30-33kd, 44-47kd
Positive IHC detected in:	Human kidney
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag4487
GenBank accession number:	BC039906
Gene ID (NCBI):	57107
Full name:	Prelyl (decaprenyl) diphosphate synthase, subunit 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.