

NSRI723WSG/0

Product Description

Ensure the quality and durability of your products with Corepoint® Scientific Temperature Test Chambers. These chambers provide a controlled temperature environment to deliver superior product testing that meets ICH Q1A standards. Improve product quality, test packaging, increase shelf life, and enhance safety with these chambers. With natural refrigerants and smart controllers, these units feature some of the most advanced technology in the industry.

Images



Certifications



General Description and Application

Storage capacity (cu. ft):	72 cu. ft.
Door:	Three glass right and left hinged, self closing, magnetic door gasket for positive seal
Shelves:	Nine stainless steel, sliding tray shelves
Drawers:	Not-applicable
Mounting and Installation:	Six swivel casters (3 1/2"), front casters locking
Interior Lighting:	Shielded, switched interior LED lighting
Airflow Management:	Plenum air distribution
External Probe Access:	Probe access port (3/4")
Insulation:	High density urethane foam cabinet and door insulation
Exterior Materials:	Powder coat steel
Access Control:	Keyed door lock
General Warranty:	Eighteen (18) months parts and labor warranty
Compressor Warranty:	Five (5) years compressor parts warranty
Product Weight (lbs):	758 lbs.
Shipping Weight (lbs):	838 lbs.
Rated Amperage:	12.4 Amps
Power Plug/Power Cord:	NEMA 5-15P conforms to UL471 requirements
Facility Electrical Requirement:	115VAC/60Hz
Agency Listing and Certification:	ETL Listed


Performance

Uniformity (Cabinet Air):	+/- 1°C from 4°C to 60°C
Stability (Cabinet Air):	+/- 0.5°C maximum from 4°C to 60°C
Temperature Control:	+/-0.1°C
Temperature Setpoint Range:	4°C to 60°C
Disclaimers:	Specifications are subject to change. For optimal uniformity at set points above 50°C, settings adjustments may be required based on load conditions and chamber configuration.

Refrigeration System	
Compressor:	Top mounted, variable speed compressor
Refrigerant:	EPA SNAP compliant, R290
RH sensor Type:	Capacitive
Temp sensor Type:	NTC
Other:	Multi-element, independent proportional electric heating, condensate evaporation

Controller, Configuration, Alarms and Monitoring	
Controller Technology:	Proportional Integral Derivative (PID) microprocessor with LCD display, °C/°F switchable
Battery Backup:	Controller
Display Technology:	LCD display
Digital Communication:	RS-485 (MODBUS)
Data Transfer:	USB port for data transfer and software updates
Chart Recorder:	Not-applicable
Adjustable Temperature Range:	4°C to 60°C
External alarm connection:	State switching remote alarm contacts
Alarms:	Visual alarm indicator, audible alarm indicator, process phase specific temp alarm, power failure alarm, sensor failure alarm
Controller Information:	Programmable ramp and soak USB data port Real time clock Event logs: 100 Process phase specific temp alarm

Dimensions					
	Width (in.)	Depth (in.)	Height (in.)	Door Swing (in.)	Total open Depth (in.)
Exterior	81"	36.25"	81.5"	22.875"	59.125"
Interior	66.75"	23.375"	59.5"		

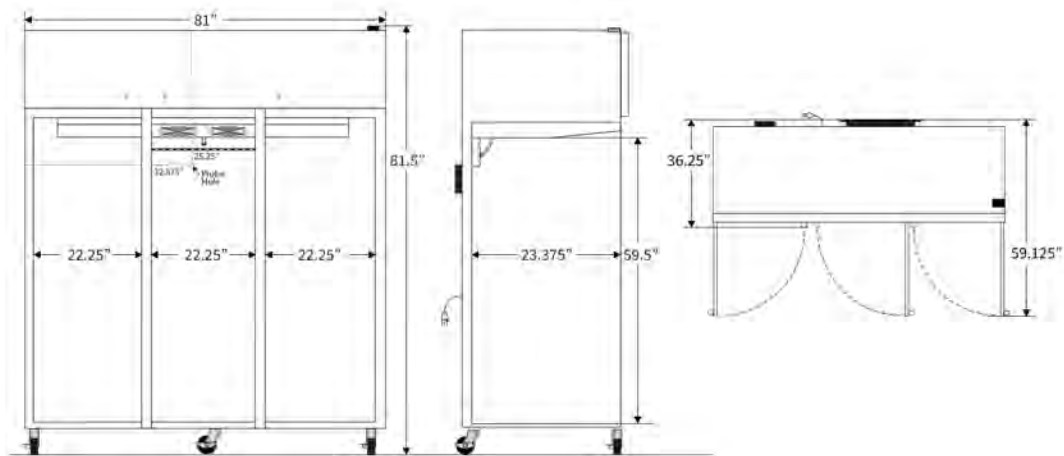

Contact		
Customer Service	800-648-4041 Option 3	customerservice@horizonscientific.com
Technical Service	800-648-4041 Option 5, Parts Option 4	technicalservice@horizonscientific.com

Rev_11252024