
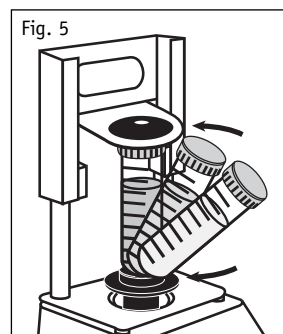
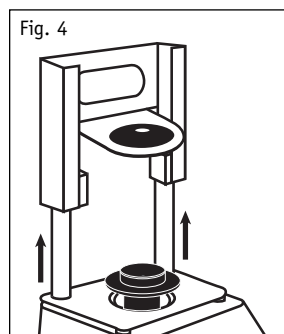
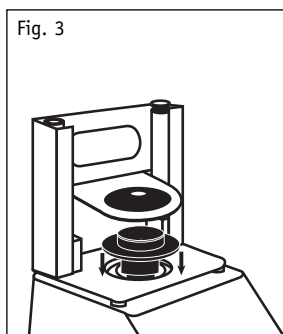
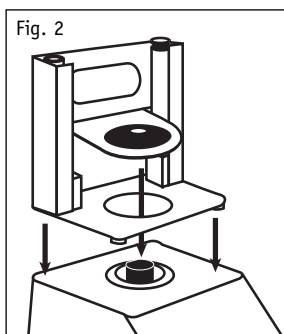


## SINGLE TUBE HOLDER INSTALLATION INSTRUCTIONS

1. Move the rocker switch to the off or standby position, in the center.
2. Remove cup head, or any other accessory from your Vortex Mixer. (Fig. 1)
3. Align the base plate of the Single Tube Holder over the post of your Vortex Mixer with the two alignment posts at the back of the unit. (Fig. 2) Lower Single Tube Holder onto Vortex Mixer, the magnets will hold the Single Tube Holder in place. To re-align this attachment, pull up from the handle until the magnets release, check the alignment and place again.
4. Now that the Single Tube Holder is in place, attach the cup head aligning the flat "alignment mark" inside the head mounting hole with the flat "alignment mark" on the post. (Fig. 3)
5. Lift up on the spring loaded upper portion of the Single Tube Holder (Fig. 4) and slide a 2.5" (6.4cm) to 4.5" (11.4cm) tube on the cup head so that the tube is aligned in the middle. (Fig. 5) The 0.5mL micro-tube holder or the 1.5-2.0mL tube holder will also fit into the Single Tube Holder. You are now ready to mix your sample.

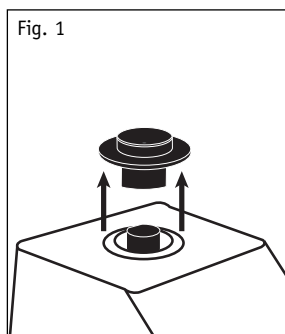

Manufactured by:

**HENRY TROEMNER LLC**

201 Wolf Drive • PO Box 87 • Thorofare, NJ 08086-0087






Phone: 856-686-1600 • Fax: 856-686-1601 • E-mail: [troemner@troemner.com](mailto:troemner@troemner.com) • [www.troemner.com](http://www.troemner.com)

715023-00

4-073-INS (REV 4/07)

## SINGLE TUBE HOLDER INSTALLATION INSTRUCTIONS

1. Move the rocker switch to the off or standby position, in the center.
2. Remove cup head, or any other accessory from your Vortex Mixer. (Fig. 1)
3. Align the base plate of the Single Tube Holder over the post of your Vortex Mixer with the two alignment posts at the back of the unit. (Fig. 2) Lower Single Tube Holder onto Vortex Mixer, the magnets will hold the Single Tube Holder in place. To re-align this attachment, pull up from the handle until the magnets release, check the alignment and place again.
4. Now that the Single Tube Holder is in place, attach the cup head aligning the flat "alignment mark" inside the head mounting hole with the flat "alignment mark" on the post. (Fig. 3)
5. Lift up on the spring loaded upper portion of the Single Tube Holder (Fig. 4) and slide a 2.5" (6.4cm) to 4.5" (11.4cm) tube on the cup head so that the tube is aligned in the middle. (Fig. 5) The 0.5mL micro-tube holder or the 1.5-2.0mL tube holder will also fit into the Single Tube Holder. You are now ready to mix your sample.

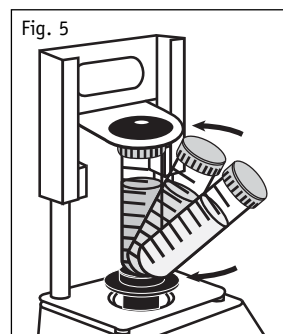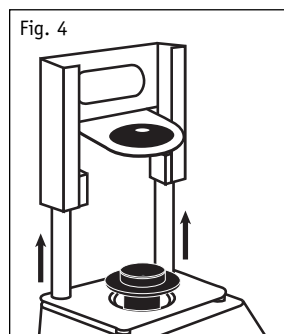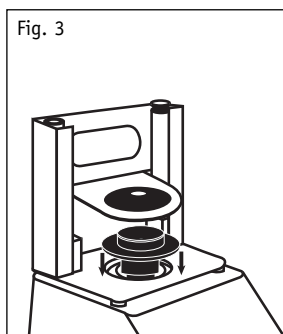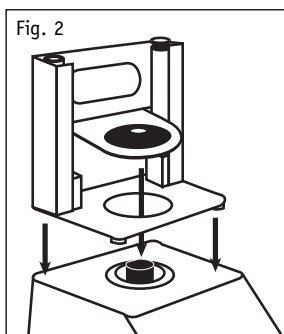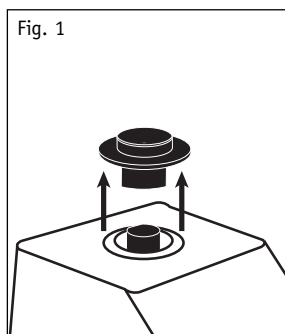

Manufactured by:

**HENRY TROEMNER LLC**

201 Wolf Drive • PO Box 87 • Thorofare, NJ 08086-0087

Phone: 856-686-1600 • Fax: 856-686-1601 • E-mail: [troemner@troemner.com](mailto:troemner@troemner.com) • [www.troemner.com](http://www.troemner.com)

715023-00

4-073-INS (REV 4/07)