

UBE2Q2 Polyclonal Antibody

Catalog number: 12581-1-AP


Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

UBE2Q2; DKFZp762C143

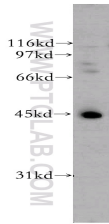

A431 cells were subjected to SDS PAGE followed by western blot with 12581-1-AP(UBE2Q2 Antibody) at dilution of 1:1000

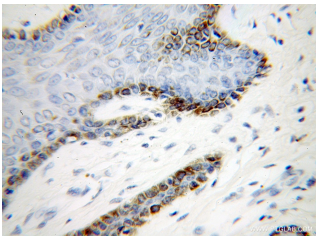

Immunohistochemistry of paraffin-embedded human skin cancer using 12581-1-AP(UBE2Q2 Antibody) at Dilution 1:100 (under 40x lens)

Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	A431 cells, mouse ovary tissue
Calculated UBE2Q2 MW:	375aa,43kd
Observed UBE2Q2 MW:	43kd
Positive IHC detected in:	Human skin cancer
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

Immunogen information

Immunogen:	Ag3276
GenBank accession number:	BC017708
Gene ID (NCBI):	92912
Full name:	Ubiquitin-conjugating enzyme E2Q family member 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.