

PLEKHH2 Polyclonal Antibody

Catalog number: 14204-1-AP

Size: 26 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

PLEKHH2; KIAA2028,

PLEKHH1L, PLEKHH2

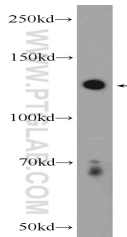

Jurkat cells were subjected to SDS PAGE followed by western blot with 14204-1-AP (PLEKHH2 Antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated PLEKHH2 MW:	168 kDa
Observed PLEKHH2 MW:	120 kDa
Positive WB detected in	Jurkat cells
Recommended dilution:	WB: 1:200-1:1000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag5409
GenBank accession number:	BC063310
Gene ID (NCBI):	130271
Full name:	Pleckstrin homology domain containing, family H (with MyTH4 domain) member 2

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.