

Mucolipin 3 Polyclonal Antibody

Catalog number: 13879-1-AP

Size: 36 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

MCOLN3; FLJ11006, FLJ36629,

MCOLN3, Mucolipin 3, TRP

ML3, TRPML3

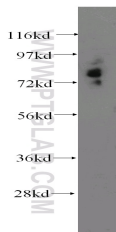

A375 cells were subjected to SDS PAGE followed by western blot with 13879-1-AP (MCOLN3 antibody) at dilution of 1:400

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated Mucolipin 3 MW:	64kd
Observed Mucolipin 3 MW:	80kd
Positive WB detected in	A375 cells
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4962
GenBank accession number:	BC060765
Gene ID (NCBI):	55283
Full name:	Mucolipin 3

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.