

DAZAP1 Polyclonal Antibody

Catalog number: 11120-1-AP


Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

DAZAP1; DAZAP1, DAZ
associated protein 1

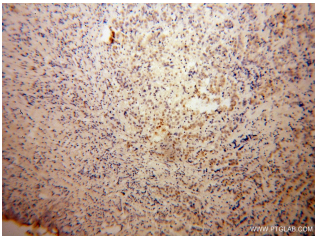

Immunohistochemical of paraffin-embedded human prostate cancer using 11120-1-AP (DAZAP1 antibody) at dilution of 1:100 (under 10x lens)

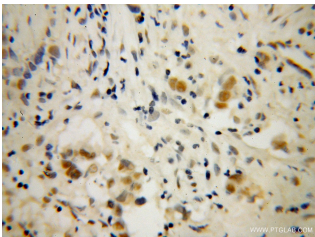

Immunohistochemical of paraffin-embedded human prostate cancer using 11120-1-AP (DAZAP1 antibody) at dilution of 1:100 (under 40x lens)

Applications

Tested applications:	ELISA, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Calculated DAZAP1 MW:	45 kDa
Observed DAZAP1 MW:	
Positive IHC detected in	Human prostate cancer
Recommended dilution:	IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1603
GenBank accession number:	BC012062
Gene ID (NCBI):	26528
Full name:	DAZ associated protein 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.