

## DPPA2 Polyclonal Antibody

Catalog number: 12689-1-AP

Size: 29 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

DPPA2; DPPA2, ECSA,

PESCRG1

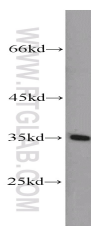

MCF7 cells were subjected to SDS PAGE followed by western blot with 12689-1-AP(ECSA antibody) at dilution of 1:500

### Applications

|                         |                                     |
|-------------------------|-------------------------------------|
| Tested applications:    | ELISA, WB                           |
| Species specificity:    | Human; other species not tested.    |
| Calculated DPPA2 MW:    | 298aa,34kd                          |
| Observed DPPA2 MW:      | 34kd                                |
| Positive WB detected in | MCF7 cells, HeLa cells, HepG2 cells |
| Recommended dilution:   | WB: 1:200-1:2000                    |

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

|                           |   |
|---------------------------|---|
| Immunogen:                | Ag3601                                  |
| GenBank accession number: | BC018070                                |
| Gene ID (NCBI):           | 151871                                  |
| Full name:                | Developmental pluripotency associated 2 |

### Product information

|                      |  |
|----------------------|--|
| Purification method: | Antigen affinity purification  |
| Storage:             | PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C. |