

PCMT1 Polyclonal Antibody

Catalog number: 10519-1-AP


Size: 28 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


PCMT1; PCMT1

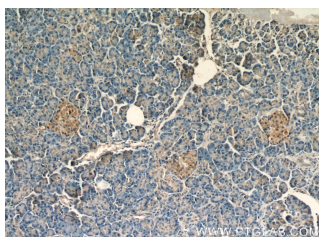

Immunohistochemical of paraffin-embedded human pancreas using 10519-1-AP(PCMT1 antibody) at dilution of 1:50 (under 10x lens)

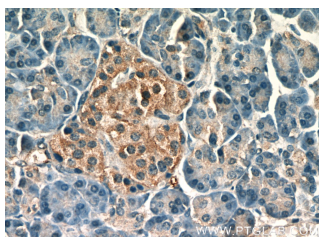

Immunohistochemical of paraffin-embedded human pancreas using 10519-1-AP(PCMT1 antibody) at dilution of 1:50 (under 40x lens)

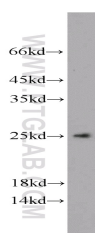

HEK-293 cells were subjected to SDS PAGE followed by western blot with 10519-1-AP(PCMT1 antibody) at dilution of 1:1000

Applications

Tested applications:	ELISA, WB, IHC
Cited applications:	WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Cited species:	Human, mouse
Calculated PCMT1 MW:	25 kDa
Observed PCMT1 MW:	25kd
Positive WB detected in	HEK-293 cells, HeLa cells, human testis tissue, Raji cells
Positive IHC detected in	Human pancreas tissue
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag0811
GenBank accession number:	BC008748
Gene ID (NCBI):	5110
Full name:	Protein-L-isoaspartate (D-aspartate) O-methyltransferase

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.