

PYGB Polyclonal Antibody

Catalog number: 12075-1-AP


Size: 28 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

PYGB; phosphorylase, glycogen;
brain, PYGB

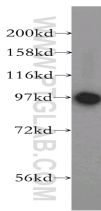

human brain tissue were subjected to SDS PAGE followed by western blot with 12075-1-AP (PYGB antibody) at dilution of 1:500

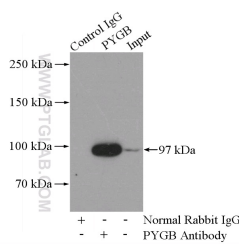

IP Result of anti-PYGB (IP:12075-1-AP, 4ug; Detection:12075-1-AP 1:500) with mouse brain tissue lysate 4000ug.

Applications

Tested applications:	ELISA, WB, IP
Cited applications:	WB
Species specificity:	Human, Mouse, Rat; other species not tested.
Cited species:	Mouse
Calculated PYGB MW:	843aa, 97 kDa
Observed PYGB MW:	97kd
Positive WB detected in	Human brain tissue, rat brain tissue
Positive IP detected in	Mouse brain tissue
Recommended dilution:	WB: 1:1000-1:10000 IP: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag2732
GenBank accession number:	BC017045
Gene ID (NCBI):	5834
Full name:	Phosphorylase, glycogen; brain

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.