

Kv4.1 Polyclonal Antibody

Catalog number: 13768-1-AP


Size: 20 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

KCND1; KCND1, Kv4.1

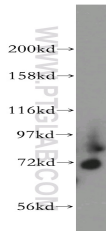

human heart tissue were subjected to SDS PAGE followed by western blot with 13768-1-AP(KCND1 antibody) at dilution of 1:300

Applications

Tested applications:	ELISA, WB
Species specificity:	Human,Mouse,Rat; other species not tested.
Calculated Kv4.1 MW:	647aa,71kd
Observed Kv4.1 MW:	71kd
Positive WB detected in	Human heart tissue, mouse brain tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag4751
GenBank accession number:	BC045659
Gene ID (NCBI):	3750
Full name:	Potassium voltage-gated channel, Shal-related subfamily, member 1

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.