

## TGS1 Polyclonal Antibody

Catalog number: 13034-1-AP

Size: 20 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

TGS1; TGS1, DKFZp762A163,

CLL associated antigen KW 2,

PIPMT, HCA137

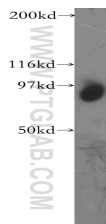

mouse testis tissue were subjected to SDS PAGE followed by western blot with 13034-1-AP (TGS1 antibody) at dilution of 1:400

### Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Calculated TGS1 MW:	852aa,90 kDa
Observed TGS1 MW:	90kd
Positive WB detected in	Mouse testis tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag3683
GenBank accession number:	BC011999
Gene ID (NCBI):	96764
Full name:	Trimethylguanosine synthase homolog ( <i>S. cerevisiae</i> )

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.