

## IL1RA Polyclonal Antibody

Catalog number: 10844-1-AP


Size: 44 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:


IL1RN; Anakinra, ICIL 1RA, IL 1ra, IL 1ra3, IL 1RN, IL1 inhibitor, IL1F3, IL1RA, IL1RN, IRAP

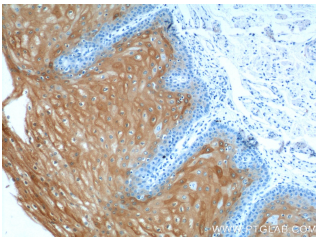

Immunohistochemistry of paraffin-embedded human oesophagus tissue slide using 10844-1-AP( IL1RN Antibody) at dilution of 1:50 (under 10x lens)

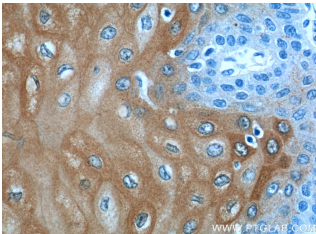

Immunohistochemistry of paraffin-embedded human oesophagus tissue slide using 10844-1-AP( IL1RN Antibody) at dilution of 1:50 (under 40x lens)

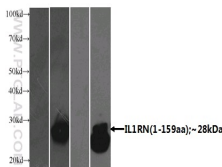

- + - + Overexpression of IL1RN-myc-his  
+ + - - Rabbit Anti human IL1RN polyclonal antibody  
- - + + Mouse Anti his-tag monoclonal antibody

Transfected HEK-293 cells were subjected to SDS PAGE followed by western blot with 10844-1-AP( IL1RN Antibody) at dilution of 1:500

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human; other species not tested.
Calculated IL1RA MW:	20 kDa
Observed IL1RA MW:	
Positive WB detected in	Transfected HEK-293 cells
Positive IHC detected in	Human oesophagus tissue, human colon cancer tissue, human stomach tissue
Recommended dilution:	WB: 1:200-1:2000 IHC: 1:20-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

### Immunogen information

Immunogen:	Ag1277
GenBank accession number:	BC009745
Gene ID (NCBI):	3557
Full name:	Interleukin 1 receptor antagonist

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.