

HNRNPA0 Polyclonal Antibody

Catalog number: 10848-1-AP


Size: 33 µg/150 µl

Source: Rabbit

Isotype: IgG

Synonyms:

HNRNPA0; hnRNP A0,
HNRNPA0, HNRPA0

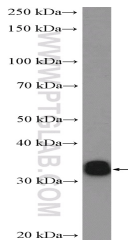

HeLa cells were subjected to SDS PAGE followed by western blot with 10848-1-AP (HNRNPA0 Antibody) at dilution of 1:600

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, Mouse; other species not tested.
Calculated HNRNPA0 MW:	31 kDa
Observed HNRNPA0 MW:	32 kDa, 34 kDa
Positive WB detected in	HeLa cells, fetal human brain tissue, human brain tissue, mouse brain tissue
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1312
GenBank accession number:	BC019271
Gene ID (NCBI):	10949
Full name:	Heterogeneous nuclear ribonucleoprotein A0

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.