

## ASB9 Polyclonal Antibody

Catalog number: 11952-1-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG

Synonyms:

ASB9; DKFZp564L0862;

FLJ20636; MGC4954

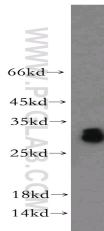

human testis tissue were subjected to SDS PAGE followed by western blot with 11952-1-AP(ASB9 Antibody) at dilution of 1:800

### Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse; other species not tested.
Positive WB detected in:	Human testis tissue, HepG2 cells
Calculated ASB9 MW:	252aa,27kd
Observed ASB9 MW:	28k-32d
Recommended dilution:	WB: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag2554
GenBank accession number:	BC001244
Gene ID (NCBI):	140462
Full name:	Ankyrin repeat and SOCS box-containing 9

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.