

NGEF Polyclonal Antibody

Catalog number: 13271-1-AP


Size: 49 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

NGEF; EPHEXIN, Ephexin 1, NGEF

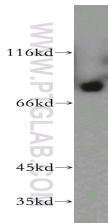

mouse brain tissue were subjected to SDS PAGE followed by western blot with 13271-1-AP(NGEF antibody) at dilution of 1:500

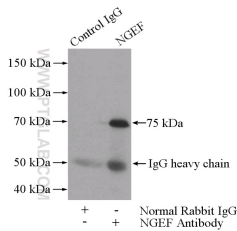

IP Result of anti-NGEF (IP:13271-1-AP, 4µg; Detection:13271-1-AP 1:500) with mouse brain tissue lysate 4000µg.

Applications

Tested applications:	ELISA, WB, IP
Species specificity:	Human, Mouse; other species not tested.
Calculated NGEF MW:	710aa,82 kDa
Observed NGEF MW:	75 kDa
Positive WB detected in	Mouse brain tissue, mouse testis tissue
Positive IP detected in	Mouse brain tissue
Recommended dilution:	WB: 1:200-1:2000 IP: 1:200-1:2000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3967
GenBank accession number:	BC031573
Gene ID (NCBI):	25791
Full name:	Neuronal guanine nucleotide exchange factor

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.