

EHD4 Polyclonal Antibody

Catalog number: 11382-2-AP


Size: 50 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

EHD4; EH domain containing 4,
EH domain containing protein
4, EHD4, HCA10, HCA11, PAST
homolog 4, PAST4

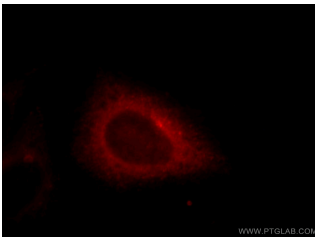

Immunofluorescent analysis of HepG2 cells, using EHD4 antibody 11382-2-AP at 1:25 dilution and Rhodamine-labeled goat anti-rabbit IgG (red).

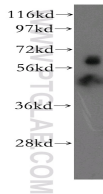

mouse testis tissue were subjected to SDS PAGE followed by western blot with 11382-2-AP(EHD4 antibody) at dilution of 1:500

Applications

Tested applications:	ELISA, WB, IF
Species specificity:	Human, Mouse, Rat; other species not tested.
Calculated EHD4 MW:	541aa, 61 kDa
Observed EHD4 MW:	61 kDa
Positive WB detected in	Mouse testis tissue, HeLa cells
Positive IF detected in	HepG2 cells
Recommended dilution:	WB: 1:200-1:2000 IF: 1:10-1:100

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag1942
GenBank accession number:	BC006287
Gene ID (NCBI):	30844
Full name:	EH-domain containing 4

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.1% sodium azide and 50% glycerol pH 7.3. Store at -20°C.