

## GALNT4 Polyclonal Antibody

Catalog number: 12897-1-AP

Size: 20 µg/150 µl


Source: Rabbit

Isotype: IgG


Synonyms:

GALNT4; GALNAC-T4;

GalNAcT4

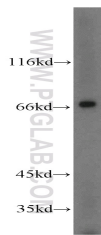

HeLa cells were subjected to SDS PAGE followed by western blot with 12897-1-AP (GALNT4 Antibody) at dilution of 1:1500

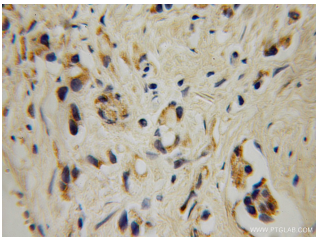

Immunohistochemistry of paraffin-embedded human prostate cancer using 12897-1-AP (GALNT4 Antibody) at Dilution 1:100 (under 10x lens)

### Applications

Tested applications:	ELISA, WB, IHC
Species specificity:	Human, mouse, rat; other species not tested.
Positive WB detected in:	HeLa cells
Calculated GALNT4 MW:	578aa,67kd
Observed GALNT4 MW:	67kd
Positive IHC detected in:	Human prostate cancer
Recommended dilution:	WB: 1:500-1:5000 IHC: 1:50-1:200

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation, FC = Flow cytometry

### Immunogen information

Immunogen:	Ag3590
GenBank accession number:	BC036390
Gene ID (NCBI):	8693
Full name:	UDP-N-acetyl-alpha-D-galactosamine:polypeptide N-acetylgalactosaminyltransferase 4 (GalNAc-T4)

### Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.