

ST3GAL6 Polyclonal Antibody

Catalog number: 13154-1-AP

Size: 28 µg/150 µl

Source: Rabbit


Isotype: IgG

Synonyms:

ST3GAL6; ST3GAL6, ST3GALVI,

SIAT10, ST3Gal VI,

Sialyltransferase 10

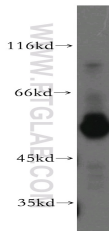

human placenta tissue were subjected to SDS PAGE followed by western blot with 13154-1-AP(ST3GAL6 antibody) at dilution of 1:200

Applications

Tested applications:	ELISA, WB
Species specificity:	Human, mouse, rat; other species not tested.
Calculated ST3GAL6 MW:	331aa,38 kDa
Observed ST3GAL6 MW:	50-55kd
Positive WB detected in	Human placenta tissue, human heart tissue
Recommended dilution:	WB: 1:500-1:5000

Application key: WB = Western blotting, IHC = Immunohistochemistry, IF = Immunofluorescence, IP = Immunoprecipitation

Immunogen information

Immunogen:	Ag3782
GenBank accession number:	BC023312
Gene ID (NCBI):	10402
Full name:	ST3 beta-galactoside alpha-2,3-sialyltransferase 6

Product information

Purification method:	Antigen affinity purification
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3. Store at -20°C.